**Rechargeable Sonic Power Toothbrush**

**Sonica B5**

**EN OPERATING MANUALS**  
KEEP THESE INSTRUCTION FOR FUTURE REFERENCE.

**IMPORTANT SAFEGUARDS**

The Power toothbrush is safe and effective when used according to directions. Like any electrical appliance, however, basic safety precautions should be followed, especially when children are present.

**DANGERS-WATER**

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not reach for a product that has fallen into water or other liquid. Unplug immediately.
- Do not use while bathing.
- Unplug the handle before cleaning and assure that it is dry before plugging it into an electrical outlet.

**DANGER-WARNINGS**

To reduce the risk of burns, electrocution, fire, or injury to persons:


- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience and knowledge, unless they have been given adequate supervision or instruction concerning use of the product by an adult responsible for their safety.
- Children should be supervised to ensure they do not play with the product.
- Use this product only for its intended use as described in this manual. This product is designed to clean your teeth, gums and tongue only. Do not use any attachments or brush heads other than those recommended by the manufacturer.
- The toothbrush is a personal care appliance and is not intended for use on multiple patients in a dental office or institution.
- If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle to avoid cracking the plastic.
- Never operate this product if it has a damaged cord, plug or charger base.
- Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- Keep the toothbrush away from heated surfaces.
- If you experience discomfort or pain in your mouth, discontinue use of this product and contact a dentist.
- Do not attempt to remove the battery.
- Do not clean the toothbrush in the dishwasher.
- The battery must be removed from the appliance before it is scrapped. The appliance must be disconnected from the supply mains when removing the battery.
- The battery is to be disposed of safely.
- Its service agent or similarly qualified persons in order to avoid a hazard.
- Regarding how to clean the tooth brush surfaces in contact with hand and skin, please refer to section <MAINTENANCE>

**MEDICAL-WARNINGS**

- If you have had oral surgery or gum surgery within the past two months, consult your dental provider before using the product.
- Consult your dental provider if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week of use.
- Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).
- Consult your physician prior to using the toothbrush if you have any medical concerns.

**SAVE THESE INSTRUCTIONS**

**CONFIGURATION**

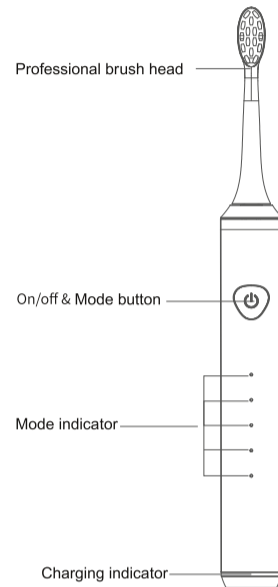**READ ALL INSTRUCTIONS BEFORE USING**

**GETTING STARTED**

**Charging the toothbrush**

Be sure to turn off the toothbrush before charging. Connect the USB cable into power supply. Indicator light will be lit, indicating the handle is charging. After the battery is fully charged, charge indicator light will turn off, indicating the toothbrush is fully charged.


**Notes:**

1. It takes at least 2 hours to charge the battery fully.
2. The toothbrush must be charged for 2 hours, if:
  - a. this is the first time using the toothbrush;
  - b. the toothbrush has not been used for more than three months;
3. During charging, the toothbrush will get warm, but this is normal.
4. The operating duration will become shorter as the battery gets older. If the toothbrush operates only for 1 short time after being fully charged, this means that the battery is nearing the end of its life cycle. Contact service center for support.

**HOW TO USE**


1. Slide the brush head onto the handle. Push brush head down until it locks into place on handle.

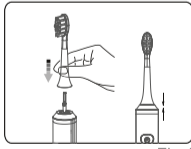Note: There is small gap between color ring and the handle.(Fig 1)

2. Wet the bristles and apply a small amount of toothpaste.

3. Place the bristles against the teeth at a slight angle towards the gum line.(Fig 2)

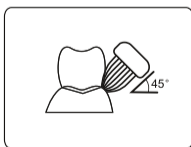4. Press the mode button to select mode. If you know the mode used last time, you need not select mode again. There are 5 modes.

- a. WHITEN: Quickly remove tooth stains and brighten your teeth.
- b. CLEAN: Effectively remove plaque and teeth superiorly cleaning.
- c. SENSITIVE: Soothing vibration for sensitive and first time users.
- d. POLISH: Polishes teeth to bring brightness naturally
- e. MASSAGE: Gently massage and protect gums

**Notes:**

5. During brushing, if press Mode button, toothbrush will work at next mode and start another 2 minutes cycle and the mode indicator will light up.
6. Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion. In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds. The brushing action will pause

between each 30 second interval, indicating when to move on to the next section. The toothbrush will automatically stop after 2 minutes of use. (Fig 3)

7. Pull brush head out from the toothbrush.(- Fig 4)
8. Rinse both the brush head and handle under running water wipe them dry. (Fig 5)

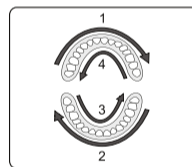
Fig 3

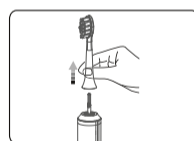
Fig 4

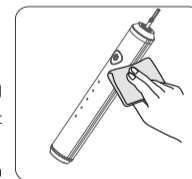
Fig 5

**MAINTENANCE**

Be sure to remove the plug from the electrical outlet before cleaning the product.

1. Remove the brush head.
2. Wipe all surface areas with a damp cloth and place the handle back into the charger. Do not submerge the handle into water.
3. Rinse the brush head and bristles after each use.
4. Do not clean in the dishwasher.

**TROUBLESHOOTING**

Problem	Possible cause	Solution
The toothbrush doesn't work.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 2 hours.
Weak or diminished vibration.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the main unit for at least 2 hours.

**DISPOSAL**





The power source of this product is Lithium-ion battery which can be recycled. At the end of the product's life, always remove the battery before you discard the product, and recycle or dispose of the battery according to local waste management requirements. If you have recycling questions, contact your local waste management office.

**Removing and recycling the battery**

Do not disassemble the product at any time other than when disposing of it. Doing so earlier will destroy the product. Observe basic safety precautions.

The appliance must be disconnected from the supply mains when removing the battery.

1. Before removing the battery, be sure to completely drain the unit of power.
2. Use a screwdriver to remove the screw cover and battery door
3. Hold the handle upside down, press the toothbrush shaft, and remove the two hooks to release the internal components of the handle.
4. Insert the screwdriver next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier.
5. Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining toothbrush components for recycling.

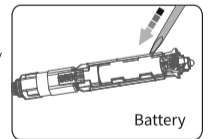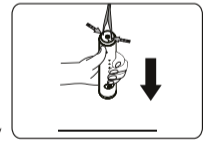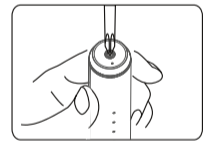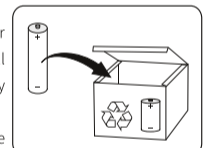Battery

Be careful not to short-circuit the negative and positive poles of the battery. Do not connect the appliance to the mains again after you have removed the battery.

Made in China for:

Global Income sp. z o.o.  
ul. Chwaszczyńska 135B  
81-571 Gdynia, Poland

www.berdson.pl

**IMPORTANT ENVIRONMENTAL INFORMATION ABOUT THIS PRODUCT**

This symbol on the device or the package indicates that disposal of the device after its lifecycle could harm the environment. Do not dispose of the unit (or battery) as unsorted municipal waste; it should be taken to a specialized company for recycling. This device should be returned to a local recycling service. Respect the local environmental rules. If in doubt, please contact your local waste disposal authorities.

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